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
USER'S GUIDE


[www.neopulse.fr](http://www.neopulse.fr)

neopulse

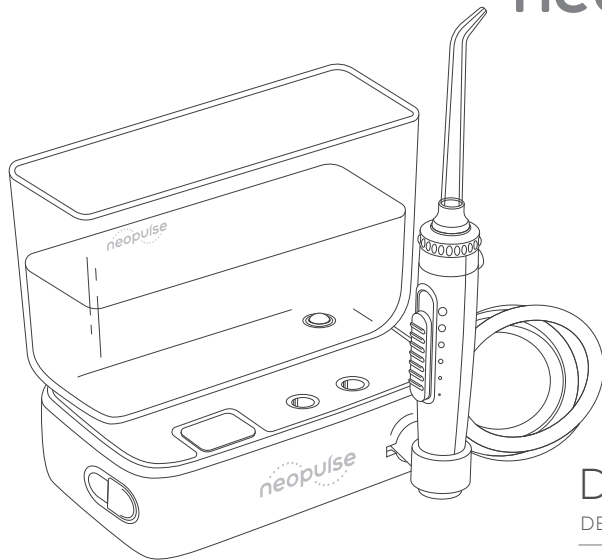JET  
DENTAIRE  
DENTAL WATER JET  
NP1 MICRO®
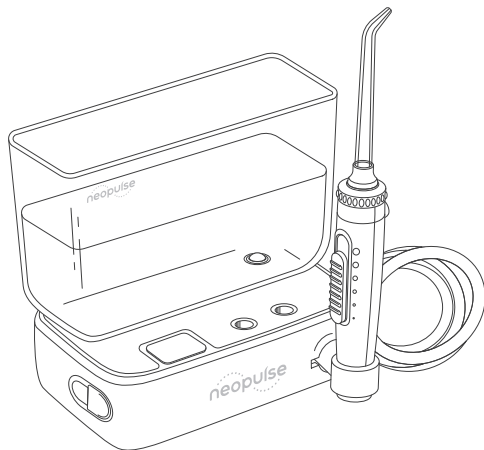
# SOMMAIRE

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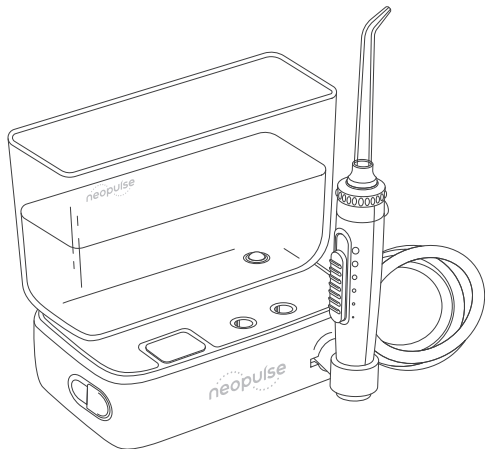
---

# PRECAUTIONS

**When using any electrical products, especially if children are present, certain basic safety precautions should always be followed, including the following:**

**BEFORE USE, READ ALL INSTRUCTIONS.**

## **DANGER**

**To avoid the risk of electric shocks:**

- Do not use this device in the bath or in the shower.
- Do not handle the USB charger with wet hands.
- Do not place the USB charger or drop into water or any other liquid.
- Check the condition of the USB charger power cord before first use, then regularly during the device life.

## **WARNING**

**To prevent against the risk of fire, burn, or injury to persons:**

- Do not connect this device to a voltage system other than the one specified on the device.
- Do not use this device if you notice that a cord or a plug is damaged.
- Do not direct the water under the tongue, inside the ear, nostrils or other delicate areas. This device is capable of producing pressure that could cause serious injury in these areas.
- Use the device as indicated in the instructions for use or follow the recommendations of your dentist, periodontist or pharmacist.
- Use only tips and accessories recommended by NEOPULSE®.
- Do not drop the device.
- Do not insert objects into the device holes.

- Keep the USB charger away from heated surfaces.
- Do not use outdoors or operate the device, where aerosol products are being used or where oxygen is being administered.
- Do not use iodine, saline solution, or concentrated insoluble essential oil in this device. The use of these products will degrade product performance and longevity. The use of these products may cancel the warranty.
- Remove any oral jewelry items before use.
- Do not use the device if you have an open wound on the tongue or in the mouth.
- If your doctor has advised you to take antibiotics before dental care, consult your dentist, periodontist or pharmacist before using this device or any other oral hygiene aid.
- Give specific instructions to children


(from 8 years) and people with special needs for proper use of this device and supervise them.

- This device does not contain any user-replaceable components, and does not require maintenance.
- Fill the tank only with lukewarm water or another solution recommended by a healthcare professional. Be careful that it does not contain dyes that could damage the system.
- If the product is opened or taken apart for any reason, the warranty will be invalidated.

#### **KEEP THESE INSTRUCTIONS.**

# PRODUCT DESCRIPTION

1. Removable tank
2. Base unit
3. Charger plug
4. Power on/off button
5. Pipe connection and handle base
6. Handle
7. Water pipe
8. Pressure control button
9. Tip ejection button
10. Directional rotary control button
11. Tank valve
12. Tip

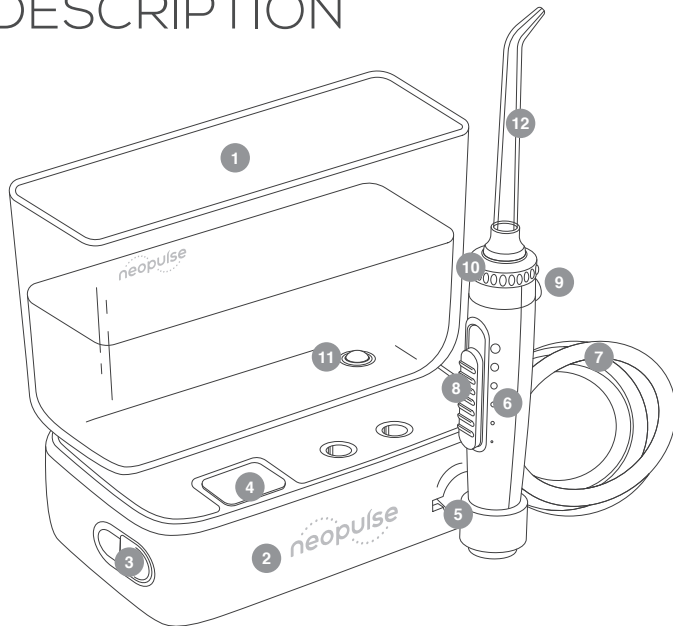TO REPLACE EVERY 6 MONTHS

---

**Classic tip**

Perfect for daily use (once or twice per day).


**Tongue-cleaning tip**

Perfect for fresh breath.

TO REPLACE EVERY 3 MONTHS

---

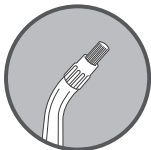**Orthodontic tip**

Perfect to clean orthodontic appliance (braces) and for general use.


**Neosoft<sup>®</sup> tip**

Perfect for deep cleaning of implants, periodontal pockets and prosthesis.


**NOTE** : replacement tips can be purchased at [www.neopulse.fr](http://www.neopulse.fr)

# FIRST USE

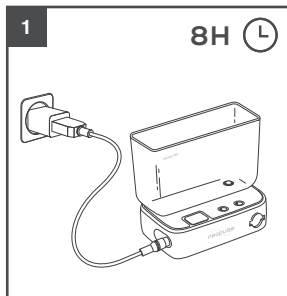
1. Charge the device before first use. Plug the USB charger into a wall outlet or into a device, which has a USB interface and place the power cord in the charger plug (3) of the device, for approximately 8 hours for the first use. The charge indicator of the USB charger changes from red to blue during the charge and to fixed blue when the device is fully charged. Do not use the device while it is in charge. Don't forget to unplug it to use it.

Under normal conditions of use, a full charge provides approximately 1 week of operation. When the battery needs to be recharged (about 5% of remaining charge), the operating speed as well as the power of the device drop. This tells the user that the device needs to be recharged. A full charge takes between 4 and 6 hours.

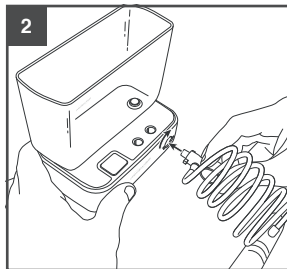
If you use the dental water jet once a day and if you are several in the family using it, it is recommended to charge it once a week.


2. Pipe / handle assembly: insert the handle base (5) into the hole in the base unit, respecting the polarizing slots. Then lock a quarter turn to the right.

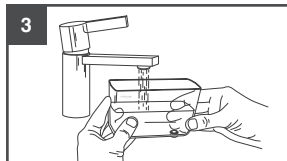
**3. Filling the tank:** remove the power cord from the device before filling the tank (1) and using the product. Remove the tank from the base unit and fill the tank with lukewarm water. Then place the tank in position on the top of the base unit and clip the valve side (11) into the base hole, press strongly. Make sure you have put the water tank in a stable position.

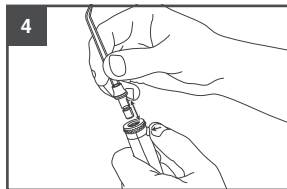**4. Tip insertion and removal:** Insert the tip (12) into the center of the rotary control button (10) on the top of the dental water jet handle. Press firmly until the tip clicks into place. To remove the tip from the handle, press the eject button (9) and remove the tip from the handle. **CAREFUL:** Do not press the tip ejection button while the device is on.

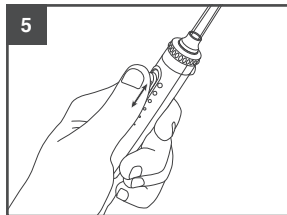
**5. Pressure settings:** Adjust the pressure by sliding the pressure control button (8), up for more pressure or down to decrease the pressure. Begin with the lowest setting and increase the pressure as preferred.

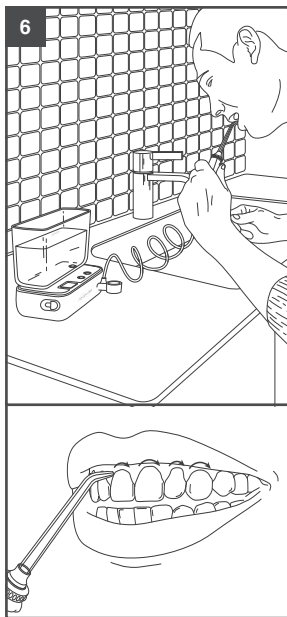
**6. Advised technique:** Lean over a sink and place the tip in your mouth. Do not place the unit too far from the sink so you do not have to pull too hard on the pipe and risk tipping over the device. Do not place the device on a surface edge to avoid the risk of it falling. Direct the tip towards the gum line (between the tooth and the gum) and switch on the appliance using the ON / OFF button (4). For best results, start with the back teeth and work towards the front teeth. Project the tip along your gum line and pause briefly between teeth inside and outside of the upper and lower teeth. Direct the water jet at a 90° angle to your gums. Close your mouth slightly to prevent water from spraying outside of your mouth, but allow water to flow freely from the mouth to the sink. When cleaning is complete, turn off the device.


**7. Use of different tips:** The Neosoft® tip is specially designed to direct water or antibacterial solutions into the furthest points of periodontal pockets (gum pockets) or under dental implants. For the use of this tip, select the lowest pressure on the device. Place the flexible end of the tip against a tooth, at a 45° angle, and gently insert the end of the tip under the gum, into the pocket, or between the implant and the gum. Turn on the dental water jet and continue moving the tip along the gum.

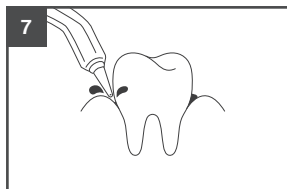**8. To use the tongue-cleaning tip**, select the lowest pressure on the dental water jet. Place the tip in the middle of the tongue, and turn on the dental water jet. Move the tongue-cleaner forwards using light pressure. Increase the pressure according to your preferences.

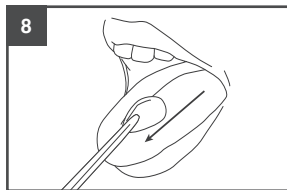**9. To use the orthodontic tip**, gently slide the tip along the gum line. Pause briefly to gently brush between the teeth and around the components of the orthodontic appliance, before moving on to the next tooth.

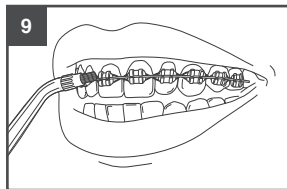**10. Use of mouthwashes and other solutions**: You can use the NEOPULSE® dental water jet for applying a mouthwash or an antibacterial solution. Following the use of special solutions, rinse the device to prevent blockages: partially fill the tank with lukewarm water and operate the device to eject all the contents of the tank into the sink (tip facing the sink).

# MAINTENANCE AND CLEANING

## MAINTENANCE INSTRUCTIONS


### After each use:

- Empty the tank completely.
- Dry the tank with a towel.
- Suck out remaining water in the device and the pipe using the tip (see diagram).


This ensures that all water is eliminated from the system in order to avoid the buildup of lime scale that will damage the unit's long-term operation.

### Maintenance and disposal:

If you need to store the dental water jet for a long time (six months or more): be sure to fully charge the unit before storing. This unit is equipped with a nickel metal hydride battery. When the battery reaches the end of its life, recycle the device to the nearest recycling center in accordance with local requirements. The battery is not replaceable. The packaging box must be placed in the paper recycling bin.

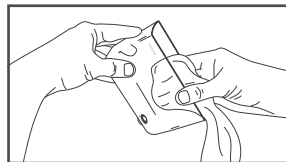Dry the tank

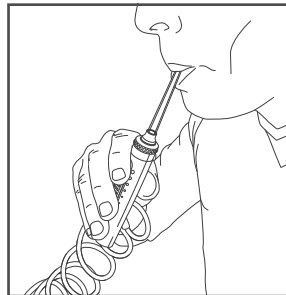Suck out the remaining water

## DEVICE CLEANING:

Clean the device when necessary. Use a soft cloth and non-abrasive cleaner to wipe the exterior. The water tank is removable, which facilitates cleaning; it can be washed in the upper basket of a dishwasher (eco program: between 45 and 55°C).

**Weekly cleaning is recommended to remove lime scale deposits.**

Deposits from hard water can accumulate in your device depending on the mineral content of your water. If the tank is not often cleaned, the lime scale may hinder the proper operation of the device.

### **Internal part cleaning:**

Add 1 tablespoon of white vinegar into a tank filled with lukewarm water. Direct the handle and the tip towards the sink. Turn the unit ON (4); let the dental water jet run until the tank is empty. For rinsing, empty the tank filled with pure lukewarm water in the same way. It is recommended that you do this once to twice a month to maximize performance.

**THIS IS A HYGIENE DEVICE,  
IT IS VERY IMPORTANT  
TO MAINTAIN IT REGULARLY  
AND KEEP IT CLEAN.**

## 2 YEAR WARRANTY

NEOPULSE® warrants that this product is free from material and manufacturing defects for a period of two years from the date of purchase.

Keep your purchase receipt as proof of the date of purchase.

NEOPULSE® reserves the right to request that the product be returned for analysis to confirm that the device is covered by the warranty.

NEOPULSE® will repair or replace the product, provided that the product has not been misused or abused, that it has not undergone any modification or damage after its purchase and that it has been used in accordance with these instructions.

For any repair or replacement service offered under the warranty, please contact Neoteam Distribution by email at [conso@neopulse.fr](mailto:conso@neopulse.fr)

Upon customer request, a copy of the warranty will be provided by email. This limited warranty does not cover accessories and consumable parts such as tips, tank, etc.

**For a better life durability of your dental water jet, NEOPULSE® is able to sell you spare parts.**

For more information, visit the website [www.neopulse.fr](http://www.neopulse.fr)  
You'll find as well video tutorials and tips.


Do not ship directly without prior agreement by email the product to the NEOTEAM DISTRIBUTION office.

Please include the model and serial numbers in all correspondence. These numbers are located under the base unit.

### Further questions?

Send us an email at  
[conso@neopulse.fr](mailto:conso@neopulse.fr)

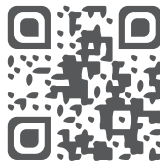[www.neopulse.fr](http://www.neopulse.fr)