THANK YOU FOR YOUR PURCHASE  
AND WELCOME TO NEOPULSE!

## SPECIAL WARNINGS

### WARNING

• NEOBABY® toothbrushes from NEOPULSE® comply with safety standards for electromagnetic devices. If you, your child, or your baby have a pacemaker or any other implanted device, consult your doctor before using this product.

• Use the device as indicated in the user manual. This toothbrush is designed exclusively for brushing the teeth and gums of babies and young children. Do not use it for other purposes.

• If your baby or child has undergone dental care in the past three months, consult your dentist, paediatric dentist, or pharmacist before using this toothbrush.

• If using the device causes discomfort or bleeding gums, consult a paediatric dentist, dentist, or general practitioner before continuing its use.

• Do not allow children to chew on the brush head.

• This device is not a toy and should not be used independently by babies and young children. Brushing must be done by an adult. Keep the device out of reach of children.

• Do not drop the device.

• Do not immerse the toothbrush in water.

• If the toothbrush gets wet, dry it thoroughly before changing the battery, and do not change batteries with wet hands. Battery replacement should only be done by an adult. Observe polarity instructions.

• Do not use a damaged battery as it may leak or explode, causing injury.

• Use a 1.5V AAA alkaline battery.

• Do not attempt to recharge non-rechargeable batteries.

• Do not expose the device to high temperatures or direct sunlight.

• This device contains no user-serviceable parts.

• If the product is disassembled for any reason, the warranty will be void.

• Do not place the brush head or handle in the dishwasher.

• Stop using the device if it is damaged.

• Replace the brush head at least every three months, or sooner if the bristles are worn. Only NEOBABY® brush heads from NEOPULSE® are compatible with the device and will not damage it.


KEEP THESE INSTRUCTIONS.

## PRODUCT DESCRIPTION

1. Baby brush head
2. Lampshade
3. Power button
4. Battery compartment
5. Suction cup

### Accessories :

- 1 baby brush head (from the first tooth)
- 1 included AAA battery (4), non-rechargeable
- 1 battery compartment cover (4)
- 1 user manual
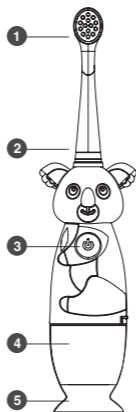
## WHAT IS THE LAMPSHADE FOR?

The lampshade (2) is used to light up and examine the child's mouth. To do this, remove the brush head (1) and use the illuminated tip (2) to spot the first emerging teeth and, most importantly, to check the quality of the brushing.

## FIRST USE:

### How to use the toothbrush?

**Step 1: Before first use**, remove the protective tab from the battery. Turn the battery compartment cover (4) counterclockwise to open it and remove the protection. Replace the battery cover by turning it clockwise.

**Step 2:** Use a fluoride toothpaste with 1000 ppm and apply an amount equivalent to a trace (for ages 0 to 2) or a rice grain (for ages 2 to 3) onto the brush head.

**Step 3:** Divide the mouth into four sections: top, bottom, left, right. This toothbrush is equipped with a 2-minute timer and vibrates every 30 seconds to signal changing to the next section.

**Step 4:** Position the toothbrush bristles across the teeth.

**Step 5:** Press the power button (3) once to turn on the lampshade (2) and use it in manual mode. Press twice to start the gentle vibrations mode. To stop the vibrations during the 2-minute brushing session, press the power button (3).

**Step 6:** Gently move the brush head over the teeth and gums using back-and-forth movements.

**Step 7:** In manual mode, change sections every 30 seconds until all four sections have been brushed. In "gentle vibrations" mode, change sections when the brush vibrates. The toothbrush automatically turns off after 2 minutes. To continue brushing, press the power button (3) again.

### How to replace the brush head?

**Step 1:** To attach the brush head, align the new head and turn it clockwise until it locks securely.

**Step 2:** Remove the brush head by turning it counterclockwise and pulling upwards.


### How to replace the battery?

**Step 1:** Hold the toothbrush upside down and turn the battery cover counterclockwise to open it.

**Step 2:** Unscrew the battery compartment with a small screwdriver.

**Step 3:** Insert a new AAA battery, observing the polarity.

**Step 4:** Replace the battery cover and tighten the screw without overtightening.

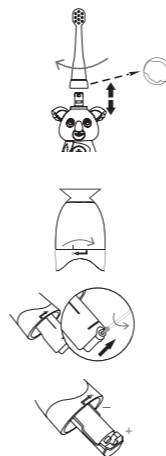**Step 5:** Align the marking line on the battery cover with the one on the brush handle.

**Step 6:** Replace the cover by turning it clockwise.

## MAINTENANCE TIPS:

**After each use**, remove any toothpaste residue:

**Step 1:** Remove the brush head and rinse the bristles with hot water.

**Step 2:** Rinse the shaft area with warm water.

**Step 3:** Clean the handle with a damp cloth.

**Step 4:** Dry the entire unit with a soft, clean towel.

### Care and disposal:

If the toothbrush needs to be stored for an extended period (six months or more), it is recommended to remove the battery. The used battery should not be disposed of in fire, nature, or with household waste; it should be collected in recycling bins designated for battery disposal. When the toothbrush reaches the end of its life, recycle it at the nearest recycling centre according to local regulations. The packaging box should be placed in the paper recycling bin. Contact us via the "support" form on our website ([www.neopulse.fr](http://www.neopulse.fr)) to give your device a second life. NEOPULSE® has a repair service (Yvelines, France) committed to extending the life of its products.

## 2-YEAR WARRANTY

NEOPULSE® guarantees that this product is free from defects in material and workmanship from the date of purchase for a period of two years, with the following exceptions. Please keep your receipt.

NEOPULSE® reserves the right to request the return of the product for analysis to confirm that the claim is covered under warranty. NEOPULSE® will repair or replace the product, provided that the product has not been subject to improper or abusive use, has not been modified or damaged after purchase, and has been used according to the instructions. For any repair or replacement service provided under warranty, please contact Neoteam Distribution by email at [conso@neopulse.fr](mailto:conso@neopulse.fr).

Upon request from the customer, a copy of the warranty will be provided via email. This 2-year warranty does not cover accessories and consumable parts such as the brush head, etc.

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Dear parents, please read all instructions carefully before the first use and keep this guide for future reference. The NEOBABY toothbrush should not be used independently by babies or young children.

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1. Remove the plastic protection from the battery.
2. To remove the toothbrush head, turn it counterclockwise, then pull it upwards.

[conso@neopulse.fr](mailto:conso@neopulse.fr)  
[www.neopulse.fr](http://www.neopulse.fr)

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