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### USER'S GUIDE

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THANK YOU FOR YOUR PURCHASE  
AND WELCOME TO NEOPULSE!

## SPECIAL WARNINGS

### DANGER

• Check the condition of the USB-C cable power cord before first use and then regularly during the life of the device. • Do not handle the USB-C cable with wet hands. • Do not immerse the USB-C cable or drop it into water. • If the toothbrush gets wet, dry it before plugging it into the power outlet to charge it.

### WARNING

• Do not connect this appliance to a system with a different voltage from that specified on the unit. • Do not use this product if the cord or plug is damaged. • Keep the USB-C cable away from hot surfaces. • NEOKIDS<sup>®</sup> toothbrushes from NEOPULSE<sup>®</sup> are comply with safety standards for electromagnetic devices. If you or your kids have a pacemaker or other implanted device, consult your physician before using this device. • Use the appliance as explained in the instructions for use, this toothbrush has been designed only for brushing the teeth and gums of children between 3 and 10 years old. Do not use it for any other purpose. • If the child has had dental treatment in the last three months, consult your dentist or pedodontists before using this toothbrush. • Give clear instructions to children on the proper use of this product, and assist with brushing up to the age of 6 - 8 years and beyond for children with special needs. • Do not drop the device. • This appliance does not contain any user-replaceable components. • If the product is opened or taken apart for any reason, the warranty will be invalidated. • Do not put the brush head, handle or charger in the dishwasher. • Replace the brush head at least every 3 months and earlier if you notice that the bristles of the brush are worn out. • Use only NEOKIDS<sup>®</sup> or NEOSONIC<sup>®</sup> brush heads from NEOPULSE<sup>®</sup>.

KEEP THESE INSTRUCTIONS.

### Step 3: Brushing

Apply the right amount of toothpaste according to the child's age (see child dental hygienist recommendations):

6 months - 3 years old: one trace

4 - 6 years old: dose of toothpaste equivalent to a pea

From 6 years old: dose of toothpaste equivalent to a chickpea

Reminder: You should help your child brush their teeth until they are 6 - 8 years old. Always check that brushing has been done correctly even if the child is independent.

Before 6 years old: The bristles of his brush are on top of his teeth.  
After 6 years old: Position the bristles on the teeth at a slight angle (45 degrees) to the teeth, along the gum line and press the On/Off and Mode button (2). For a uniform brushing, divide the child's mouth into 4 separate zones and spend 30 seconds brushing time per zone (4 zones X 30 seconds = 2 minutes).

With Smart Brushing, every 30 seconds, a gentle vibration will tell you to change the brushing zone.

The NEOKIDS<sup>®</sup> toothbrush stops automatically after 2 minutes of brushing cycle, thus complying with the recommendations of dentists and pedodontists.

Recommendations for sonic electric brushing:

To clean the outer surfaces: brush top and bottom separately. Move the brush in a rotating motion, always from pink to white, to brush your teeth and gums without exerting pressure so as not to damage the bristles too quickly.

To clean the inner surface of the teeth: make several vertical brushings on each tooth, always from the gum to the tooth.

Continue for the duration of the timer (2 minutes), until the toothbrush stops.

When using a sonic toothbrush for the first time, it is possible to at the beginning to feel some tickling. This feeling is common and disappears quickly after several uses.

Our extra advice:

Once the 2-minute brushing cycle is complete, you can go back and insist on areas prone to stains and cavities, but also on the chewing surfaces of the teeth and in crannies.

Change the child's brush head at least every 3 months.

NEOKIDS<sup>®</sup> toothbrushes are equipped with ultra-soft brush heads. Their rounded nylon bristles, are designed to effectively remove plaque and take care of your child's oral sensitivity. The colored bristles in the center serve as a visual guide for toothpaste dosing.

Important: make sure that the child does not chew the brush head.


## PRODUCT DESCRIPTION

1. Kid's brush head
2. On / Off
2. 2 brushing modes
3. Brushing mode indicators, 2 modes:
  - white light (soft brushing)
  - blue light (extra-soft brushing)
3. Battery level indicator
4. USB-C charger plug
5. Suction base

### Accessories:

- 1 small ultra-soft brush head (0,127 mm)
- 1 USB-C cable
- 1 traveling pouch
- A set of stickers for the child to personalize his toothbrush

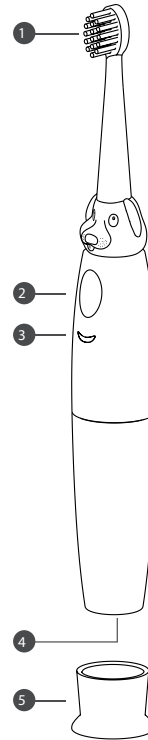All NEOPULSE brush heads (children and adults) fit NEOKIDS<sup>®</sup> and NEOSONIC<sup>®</sup> toothbrush handles. (for children, use the sensitive or gum care mode for the NEOSONIC<sup>®</sup> handle).

From the age of 6 years old you can start educating your child about interdental cleaning (cleaning the spaces between the teeth) by completing his brushing once a day with one of our NEOPULSE<sup>®</sup> dental water jets, in low power (soft), because a tooth has not only 3 sides but 5! Completing the child's brushing with a dental water jet helps to limit the risk of interdental cavities.

As soon as 2 teeth touch each other, complete the brushing every evening with NEOPICK<sup>®</sup>, an easy-to-use dental floss thanks to his handle.

The NEOKIDS<sup>®</sup> toothbrush is safe for orthodontic appliances. The brush head is likely to wear out more quickly, it is recommended to change it every 2 months.

## MAINTENANCE ADVICE

After each use, remove all toothpaste residues:

Step 1: remove the brush head and rinse the bristles well with warm water

Step 2: rinse the metal shaft area with warm water

Step 3: clean the handle with a damp cloth

Step 4: dry the whole area with a soft, clean towel

### Maintenance and disposal:

If the child's toothbrush is to be stored for an extended period of time (six months or more), be sure to fully charge it before storing it in a cool, dry place. This device is equipped with a lithium-ion battery.

When the battery reaches the end of its life, recycle the toothbrush at the nearest recycling centre in accordance with local requirements.

The packaging box must be placed in the paper recycling bin. You can also contact us via the 'assistance' form on our website ([www.neopulse.fr](http://www.neopulse.fr)) so that we can give your device a second life. NEOPULSE has a repair department (Welnes department in France) dedicated to extending the life of its products.

## 2 YEAR WARRANTY

NEOPULSE<sup>®</sup> warrants that this product is free from material and manufacturing defects for a period of two years from the date of purchase. Keep your purchase receipt as proof of the date of purchase.

NEOPULSE<sup>®</sup> reserves the right to request that the product be returned for analysis to confirm that the device is covered by the warranty. NEOPULSE<sup>®</sup> will repair or replace the product, provided that the product has not been misused or abused, that it has not undergone any modification or damage after its purchase and that it has been used in accordance with these instructions. For any repair or replacement service offered under the warranty, please contact Neotem Distribution by email at [conso@neopulse.fr](mailto:conso@neopulse.fr)

Upon customer request, a copy of the warranty will be provided by email. This limited warranty does not cover accessories and consumable parts such as brush heads, USB-C cable, etc.

## FIRST USE:

### How to use the toothbrush?

1. Mandatory, charging the device before first use


Connect the USB-C cable to a universal power adapter and place the power cord to the charging plug (4) of the device, for 6 hours for the first use (even if your toothbrush arrives already charged). The battery level indicator (3) on the device changes from flashing green to a fixed green when the toothbrush is fully charged.

Do not operate the toothbrush while it is in charge. Don't forget to unplug it to use it. Under normal conditions of use (1 child brushing twice a day for 2 minutes), a full charge allows 1 month of use. When the battery needs to be recharged, the battery level indicator (3) will flash red. A full charge takes about 3 hours.

### 2. Use of the toothbrush:

#### Step 1: Brush head set up

Align the brush head so that the bristles are in line with the front of the handle, then slide it along the shaft and clip it in place. It is normal for there to be a small gap between the brush head and the handle, the gap ensures that the toothbrush does not stick to the motor during high frequency vibration and thus allow the motor to be at the right power.

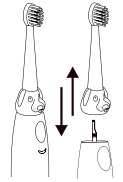To remove it: pull in the axis of the brush and do not shear the brush, as this can damage the brush motor (the warranty will be invalid).

#### Step 2: Customise the toothbrush with the child and choose from the 2 modes.

Press the setting button (2) and select the brushing mode (3), that means the vibrations amplitude that best suits the child's teeth and gums:

- White light: gentle brushing
- Blue light: extra gentle brushing

Always start with the softest mode and then find the one that best suits the child's need. Once the toothbrush has been customised, the chosen setting is stored in memory. Simply press the ON button (2) once to start brushing.

To stop the toothbrush, press the setting button (2) again for 5 seconds.

#### Last Step: Toothbrush customization

Add the finishing touch with the child by sticking the stickers they like on their toothbrush.

neopulse

[www.neopulse.fr](http://www.neopulse.fr)